

M A X I C O M F O R T TM S E R I E S



Rise







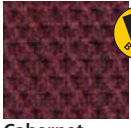

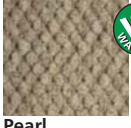



Available NOW!

Featuring Pocket Coil Seat Springs!



Cirrus Medium
PR-508M

MAXI COMFORTTM
Ultimate Recline Technology

-  **Palomino** 
-  **Admiral** 
-  **Evergreen** 
-  **Cabernet** 
-  **Pearl** 
-  **Brisa Buckskin*** 

* upcharge applies

*Cirrus*TM
by Golden

The CirrusTM with MaxiComfort lift chair has taken comfort to a whole new level! Our pocket coil spring technology and supportive foam casing molds to your lower body, gently cradling you in plush comfort. Experience the newest heights in lift chair comfort – experience the Cirrus!



Rise. Roll. Ride.

GOLDEN
technologies

PROUDLY

MADE IN AMERICA
with Domestic & Imported Materials

YOUR LIFE IN MOTION

Model Number	PR-508M
Size	Medium
Weight Capacity	375 LB
Recline Positions	Maxi Comfort
Standard Fabrics ¹	6
Fabric Options	All
Chaise Seat	Yes
Coil Springs	Yes
Back Type	Biscuit
Head Rest Cover	Yes
Arm Covers	No
Battery Back Up	Yes
Overall Width	32"
Overall Height	43"
Width Between Arms	20"
Floor to Top of Seat	20"
Seat Depth	20.5"
Seat to Top of Back	29"
Seat to Top of Arm	5.5"
Seat Angle Raised	15°
Floor to Top of Foot Rest Reclined	23"
Floor to Top of Back Reclined	22"
Back Angle Closed	100°
Back Angle Reclined	170°
Foot Rest Angle Closed	90°
Foot Rest Angle Opened	180°
Distance Required from Wall Reclined	26"
Extension Reclined from Back of Seat to End of Foot Rest	42"
Life Time Warranty, Chair Frame, Lift Frame, Recline Mechanism	Yes
Warranty on Electrical parts	3 YR
Pro Rated Years 4 Thru 7	Yes
Warranty on Mechanical Labor	3 YR
Warranty on Electrical Labor	1 YR
Removable Back Option Available	No
Smart Tek System Standard	Yes
Customizable *	No
Left Hand Control Option Available	Yes
Extra Pocket Option Available	Yes
Heat & Massage Option Available	Yes
Gel Memory Foam Option Available	Yes
Foot Rest Extension Option Available	Yes
Straight Lift Option Available	No
User Height (suggested)	5'4" -5'10"

¹ Buckskin Brisa is a standard, stock fabric on PR508M. Upcharge applies on all sizes.

* See Customize Chart for Dimension limits
 Sizes may vary depending on fabric, filling material, upholstery and carpet thickness. All measurements taken on concrete floor using levels and metal rulers.

Hand-blown, Overstuffed, Conjugated Microfiber Dacron Fill
SUPER SUPPORTIVE

Quick Release AutoDrive Hand Control
EASY DISCONNECT

Widest Fabric Selection in Industry
FITS IN ANY HOME DÉCOR

Adjustable Anti-Skid/Scrape Leg Levelers
STABILITY & SUPPORT

Kiln-dried Northern Maple Hardwood
 • Double-doweled at Stress Points
 • Reinforced Miter Blocks
 • Stapled and Glued
STRONGEST IN THE INDUSTRY

10 Gauge Steel Scissor Mechanism
STRONG & LONG LASTING

Rear Landing Legs
EXTRA STABILITY

5 Heavy Duty Springs
LONG LASTING FOR INCREDIBLE COMFORT



Our patented technologies, combined with our superior kiln-dried northern maple hardwood frame design, provide the most comfortable lift and recline experience.



The Cirrus™ features a combination of pocket coil seat springs surrounded by foam that molds to your lower body for the ultimate support and comfort.



Battery Backup System

Batteries not included. Battery Backup System is designed to bring the chair from the reclined to the seated position in the event of a power outage.



Smart Tek

Our luxurious lift and recline chairs are equipped with electronics to help assure a secure lift every time you need it.



Literature is current at time of printing. Golden Technologies reserves the right to make changes to the product or literature at any time.



Complies with UL 73 and CSA C22.2 No. 68
 Listing # E112417
 FNAE 80-214A

